ABS MASTER PRO SUPER

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- **Thank you for purchasing our fitness equipment**
- **This product is manufactured under strict control. You can use it safely**
- **To avoid damage or loss caused by improper use of the equipment, please read the "Precautions" section. Before you start training, make sure you fully understand how to use it**
- **Please keep the manual**

PRECAUTIONS

Before use, read the precautions and use the equipment correctly and safely. Precautions in the manual are intended to prevent hazards, damage to the property or the user and other people.

- 1 Allergists
- 2 Skin problems
- 3 Blood pressure problems
- 4 (Users) taking drugs
- 5 Patients who have a malignant tumor
- 6 Sensory brain disorders (e.g. epilepsy, etc.).
- 7 Sensory disturbances caused by diabetes
- 8 Patients with chronic pain such as back pain, sprains, pulled muscles, etc.
- 9 Patients with osteoporosis
- 10 Patients with fever exceeding 38C
- 11 Patients requiring increased rest
- 12 Patients who have skin problems after using gel sheets
- 13 People not feeling well
- 14 People waiting for a medical procedure

WHEN USING GEL SHEETS PLEASE OBSERVE THE FOLLOWING SYMPTOMS, OTHERWISE IT MAY CAUSE PAIN OR BURNING:

- Do not use damaged, degraded or dirty gel sheets
- Do not bend or fold the gel sheets as this may damage them
- Please place the gel sheets precisely on the skin. If the fit with the skin is not accurate the pulse does not flow evenly through the sheet. Consequently, it may cause stronger stimulation in one place
- Please turn off the device before removing the sheets
- Make sure to connect the gel sheets to the electrode. Do not touch the electrode directly while charging.

DO NOT USE THE DEVICE IF YOU BELONG TO THE GROUP OF PEOPLE LISTED BELOW. USING THE DEVICE MAY CAUSE AN ACCIDENT AS WELL AS PHYSICAL OR SKIN PROBLEMS:

- People with heart disease
- Pregnant, postpartum women
- People with heart assist devices such as a pacemaker
- People with life support devices
- Persons using electronic devices for medical examinations, e.g. holter


PLEASE OBSERVE THE FOLLOWING SAFETY PRECAUTIONS WHEN USING THE EQUIPMENT, OTHERWISE AN ACCIDENT, PERSONAL INJURY, ELECTRICAL SHOCK, DAMAGE OR SOILING OF THE PRODUCT MAY TAKE PLACE

- Do not place the device or gel sheets in a microwave or near an ignition source
- Do not use the device for purposes other than those described in the manual
- Do not apply on children or animals
- Do not consider it as a toy for children or pets
- Don't spill liquid on the equipment
- Do not use during exercise or while driving
- Do not use in a warm room such as a bathroom

DO NOT USE THE PRODUCT ON THE FOLLOWING BODY PARTS, OTHERWISE IT MAY CAUSE AN ACCIDENT AS WELL AS SKIN OR PHYSICAL PROBLEMS

- The peripheral part of the chest and shoulder blades
- Above the neck (head, face)
- Do not use on wounds and mucous membranes, eyes, mouth and genitals
- In areas, where hair removal procedure has been performed, do not use on the day of the procedure
- Where the plastic surgery was performed
- In a place with implants
- Spine
- In a place with a tattoo
- On the abdomen during the menstrual cycle

PRODUCT SPECIFICATION

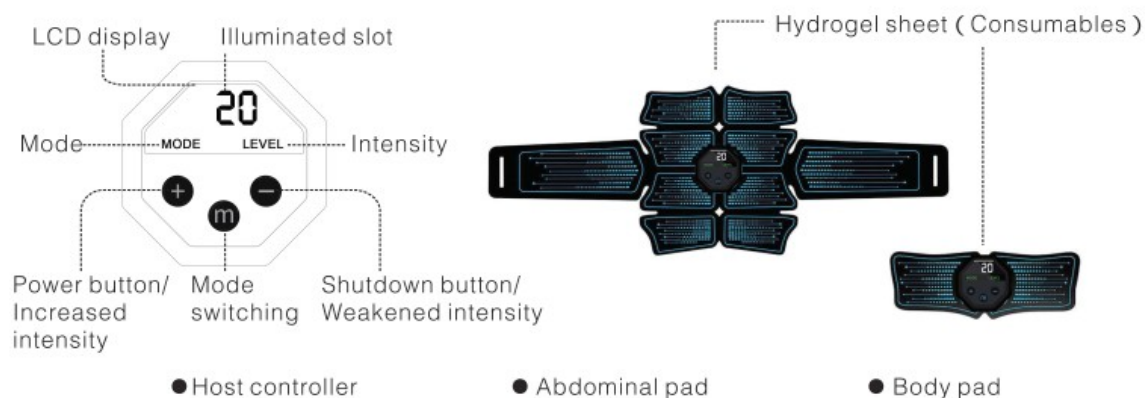Product name		ABS MASTER PRO SUPER ELEKTROSTYMULATOR HMS	
Product number	ABS	Dimensions	50*48*12MM
Weight	17g	Power	DC3.7V, 90mAh
Output intensity	0-20 levels	Levels	1-10
Materials	ABS, PC	Tested	CE
Operating time	20 min., automatic shutdown after 1 min of rest		
Recommended ambient temperature	10-40C, Ambient humidity: 30-85%		
Accessories:	Remote control/controller, gel sheets, USB cable, manual		

	Gel sheet for the abdominal pad	Gel sheet for the shoulder pad
Weight	1 pc. 10g	2 pc. 20g
Size	190(l)x160(w)x2,2(h)mm	205(l)x77(w)x2.2(h) mm
Materials	ABS, PC	ABS, PC
Quality Inspection	CE	CE


In order to improve the product some specifications may be amended, we are sorry for any inconveniences.

The gel sheet is a wearable item. If the sheet is clearly worn or damaged, a new one should be purchased.

Charger settings

- Connect the USB cable to the adapter then to the power supply; total charging time equals 2 hours.
- Red light informs on active process, light turned off indicates that charging has been completed
- The product does not include an adapter / socket tip, but can be plugged into any standard USV adapter (Output power: 5V)

all home USB power adapters (Voltage output : 5V) .

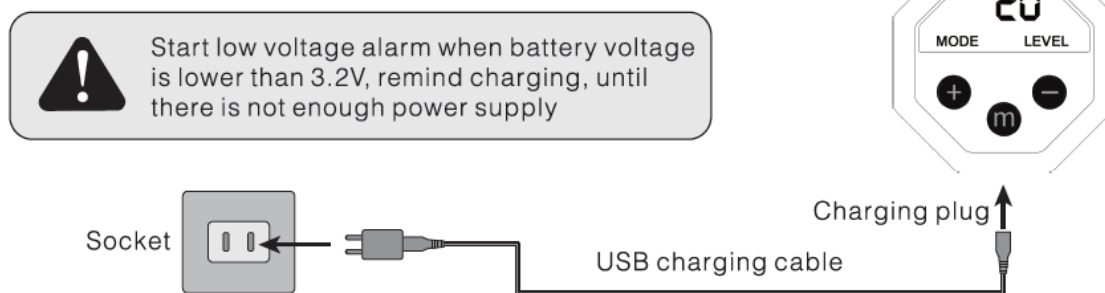HOW TO USE GEL SHEETS

- The abdominal and shoulder sheets are installed in the same way
 - Below find the correct way to apply the pad, the adhesion and the length of the device can be shortened due to excessive wear and tear
1. Connect the controller to the abdominal pad. Align the electrode hole with the electrode port
A sheet is relatively soft, be careful and prevent the two ends of the sheet from sticking together
 2. Remove the protective film from the pad
 3. Place the pad on the body
 4. After a workout has been completed put the pad on the smooth side of the transparent film to protect the gel sheets

Basic manual

Recommendations for the amount and intensity of training: It is recommended to train once a day. Please reduce the intensity when doing the first workout. Increase the intensity gradually after a slow adaptation

Place the pad on the body where the training is to be addressed to


Key principles:

- Remove sweat from the body before applying the pad
- Try not to touch the gel sheets
- Fit the gel sheets to your body. If there is a gap between the gel sheet and your body, the stream shall not flow evenly which may cause local tingling and pain

Abdominal pad - place the gel sheets as in the diagram below. If the stream is too weak, move the pad up or down by about 1-2 cm. Turn off the device before changing placement of the pad.


Shoulder pads - waist, arms, legs

1. Turn on the device

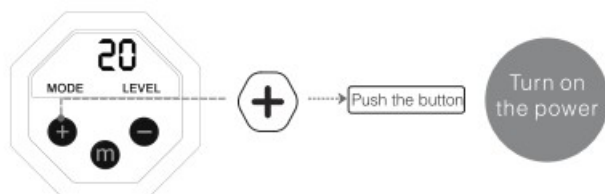If after turning on the device you find that the gel sheet does not fit on your body accurately, you may experience local tingling or pain


2. Adjusting the intensity of your training

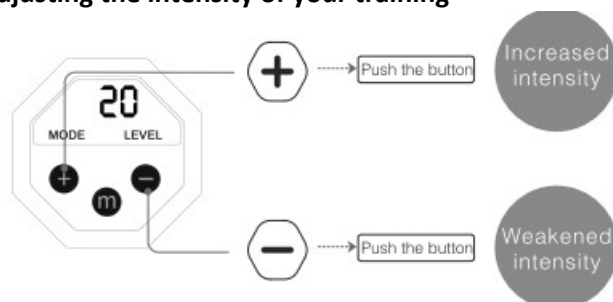- When the device is turned on press the + button, the blue light will flash and indicate the intensity of the workout at the level 1
- The training intensity strength may be selected from 20 levels (1-20)
- During the first use start with the lowest intensity level and gradually increase the intensity after a previous adaptation to a lower level
- Power of the device is consumed during a training which may translate into a weaker stream sensation. Be sure to charge the battery.

3. Program change

4 Mode conversion

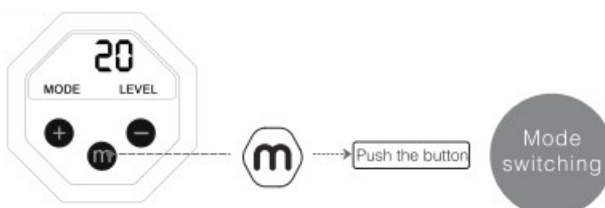- The training is divided into 10 types (1 - 10)
- Each type automatically shuts off after 20 minute operation
- Re-set to 20 minutes after switching modes starts intensity from the zero (level)

Mode1
Mixed training

Mode2
Fat warming up

Mode3
Muscle warming up

Mode4
Jogging

Mode5
Swim

Mode6
Push ups

Mode7
Sit-ups

Mode8
Barbell

Mode9
**Fat accelerated
burning**

Mode10
**Muscle intensive
training**

- Mode 1 - Mixed Training**
- Mode 2 - Fat Burning - warm-up.**
- Mode 3 - Muscle warm-up**
- Mode 4 - Jogging**
- Mode 5 - Swimming**
- Mode 6 – Push-ups**
- Mode 7 - Squats**
- Mode 8 - Barbell**
- Mode 9 - Intense fat burning**
- Mode 10 - Intensive muscle building**


Shutting down the device

- The device runs for 20 minutes; then it will automatically shut down; the blue light goes out
- The device will automatically turn off if there is no operation for one minute in the stand-by mode
- To turn off the device abruptly during a workout, hold down the off button (" - " button) for 1 second. The device shuts off.

3. Removing (taking off) the device from the body after training


Proper removal of the device:

Grasp the end of the pad and gently pull it from top to bottom as shown in the picture below:

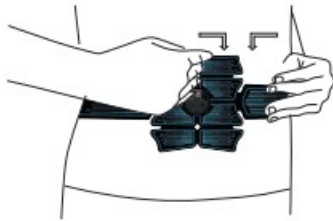Incorrect device removal:

Don't pull the pad too hard, don't try to grab small parts to pull the belt off. You may damage the gel sheets.

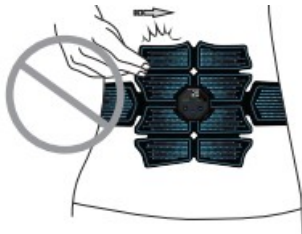4. Put the belt on the protective film to protect the gel sheets. Put the device into your bag, store safely until your next workout

MAINTENANCE

5. To protect from scratches, gently wipe away sweat and moisture with a dry cloth
6. If the device is dirty, soak a soft cloth in a warm water with a neutral detergent. After wringing it out, gently wipe the dirt from the surface
7. Soak a soft cloth in clean water, then wring and wipe the detergent from the pad
8. After cleaning, dry the device in a natural, shaded well-ventilated place

REPLACEMENT OF THE GEL SHEETS

Please follow the three steps below:

4. Remove the blue film from the gel (the blue side is the front)
5. Place on a silver electrode
6. After carefully placing the sheet on the pad remove the transparent film from the back of the gel sheet (the transparent film is the back)

Note: Remove the film gently, too much force may damage the gel sheet.

Do not use if the gel sheet is damaged, dirty or its condition has deteriorated.

Solutions

Status:	Probable cause	Solution
No stimulation	Device turned off	Press mode button "+"
	Gel sheet does not adhere to the skin	Match the sheet directly to the skin
	Gel sheet is damaged	Replace the sheet
	The gel sheet is dirty, dusty	Replace the sheet
Stimulation is too low	Low intensity level	Increase the intensity with the + button
	The power of the device is low	Charge the device
	Damaged gel sheet	Replace the gel sheet
	The gel sheet is dirty	Replace the gel sheet
Stimulation is too strong	High level setting	Reduce training intensity
	Sheet does not adhere to the skin	Make sure the sheet adheres to the skin
Sudden obstacle to the stimulation	Remaining power is exhausted	Charge the device
Gel sheet is not sticky	Gel sheet is dirty	Replace the gel sheet
	Gel sheet is damaged	Replace the gel sheet
Redness and pinching	The sheet does not adhere to the skin	Make sure the sheet adheres to the skin
	Gel sheet is dirty or dry	Replace the gel sheet
	Gel sheet is damaged	Replace the gel sheet

If, after these solutions have been applied and the problem still exists, contact your product vendor.

Even though we do our best to ensure the highest quality of our products, single defects or omissions might occur. If you note a defect or missing part, we invite you to contact us.

IMPOERTER: **ABISAL SP.Z.O.O.**,

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abisal@abisal.pl

www.abisal.pl

Marking the equipment with a symbol of a crossed-out waste container informs about prohibition of placing used electrical and electronic equipment with other household waste. According to the WEEE Directive on used electrical and electronic waste management, separate disposal methods must be applied for this type of equipment.

A user, who intends to dispose this product, is required to deliver it to a collection point for used electrical and electronic equipment in order to contribute to reuse, recycling, or recovery and consequently, to protect the environment. To do this, contact the point where the device was purchased or local government officials. Hazardous components contained in the electronic equipment can cause long-lasting adverse changes in the environment, as well as have a harmful effect on human health.

GUARANTEE CARD

Article name:.....

EAN code:.....

Date of sale:

GUARANTEE TERMS:

- The Seller provides guarantee on behalf of the Guarantor within the territory of the Republic of Poland for the period of 24 months from the date of sale: electric power assisted booster 24 months, charger 6 months.
- The Guarantee will be recognised by the shop or service centre after the client provides:
 - clearly and correctly filled-in guarantee card with the sale stamp and the seller's signature
 - valid purchase confirmation for the equipment including the date of sale / bill, complained product.
- Possible defects and damages revealed during the guarantee period shall be remedied free of charge within not more than 21 days from the delivery of the product to the shop or service.
- Should it be necessary to import some components, the guarantee term can be extended by the period required to import such part, however such period shall not be longer than 90 days.
- The guarantee does not cover:
 - mechanical damages and subsequent defects,
 - damages and defects resulting from incorrect storage and usage of the equipment against its purpose,
 - incorrect assembly and maintenance,
 - damages and consumption of such elements as: cables, straps, rubber parts, pedals, sponge grips, wheels, bearings, etc.
- The guarantee expires in the event of:
 - expiry of the validity term,
 - self-repairs,
 - failure to follow the rules of correct operation.
- The product handed over for repair should be complete and clean. If some parts are missing, the service is entitled to refuse to accept the product for repair. The service may refuse to accept a dirty product or clean in on the client's expense after the client's written permit.
- The guarantee does not cover activities related to assembly, maintenance, which, according to the instruction, the user is obliged to perform themselves.
- The Guarantor also informs that they offer post-guarantee service.
- The product should be provided in its original packaging and should be secured for shipment.
- To exercise the guarantee follow the procedure posted on the website: <https://serwis.abisal.pl/>

In case of non-conformity of the sold thing with the contract, the buyer is entitled by law to legal remedies from and at the expense of the seller. The guarantee does not affect such remedies.

THE EQUIPMENT IS NOT INTENDED FOR REHABILITATION AND THERAPY

NOTES ON THE COURSE OF REPAIRS

Item	Date of notification	Date of provision	Course of repairs	Signature of the recipient (shop, owner)

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